

“Slide labeling alone accounted for 67% of defects”

Source: The Henry Ford Production System: Measures of Process Defects and Waste in Surgical Pathology as a Basis for Quality Improvement Initiatives ©2007 American Society for Clinical Pathology

# YOUR NEW LAB BUDDY

# pslim

Laboratory Slide Printer  
by  
**Accuplace**



*Compact Footprint*



*Works With Any Slide Or LIS*

*Industry Leading Speed*

*No Computer Needed*

*Most Cost Effective*

accuplace.com | (954) 791-1500

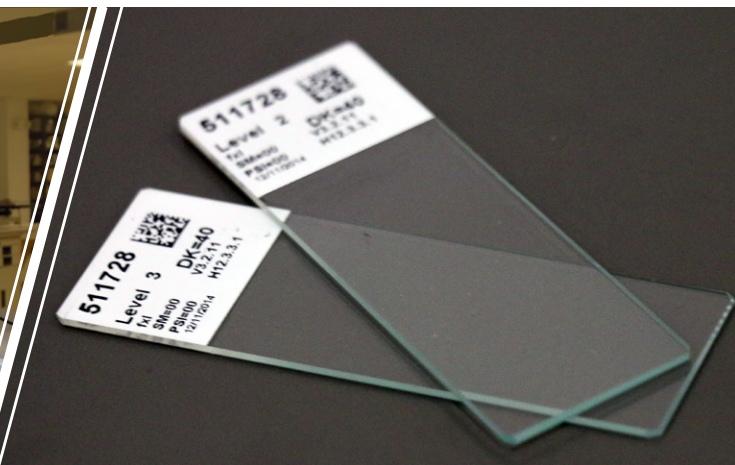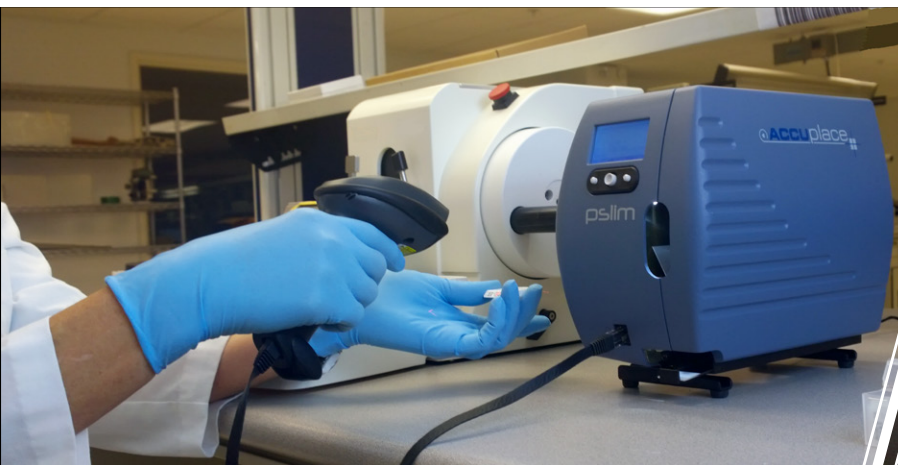

## Specifications

Printing Registration Accuracy	± 0.002"/0.05 mm*
Size:	
Width	6.0"/152 mm
Height	8.9"/226mm
Depth	10.4"/264mm
Weight	9.5 lbs/4.3 kg
Input Voltage	100-240 VZC 50-60Hz 1ø
Power Requirements	Universal 100-240 VAC 47-63Hz 1.2A max.
Color	Grey/Blue
Slide Type	3"x1" or 25mm x 75mm x 1mm w. frosted or color end*
Media Type	Thermal transfer ribbon, 1" wide roll
Media Capacity	1 roll, approximately 10,000 Slides
Media Color	Black
Throughput	3-12 Slides/minute typical, slide dependent

\*Slide dependent.

## Additional Features

Recessed connector area allows machine to go against a counter back-splash or even in a corner

Display shows data ready for printing

User can select number of slides to print

Machine can automatically print slide number/number in series if desired. (1 of 4, 2 of 4, etc.)

Format controls allow automatic printing of fixed information such as lab name, username, etc.

Highly ergonomic design suitable for counter level or shelf level operation

Slides are held inside the unit's side door to protect them from dust and the environment

UL/CE marked universal AC power adapter included


Ethernet and USB interface for PC and LIS integration

Powerful microprocessor enhances the machine's speed and ease of setup

Multiple fonts can be printed on each slide - large & small readable characters and barcodes

Front data connector interfaces with a barcode scanner or a keyboard/scanner adapter

Multi-language menu system (English, French, Italian, German, Chinese & Japanese)

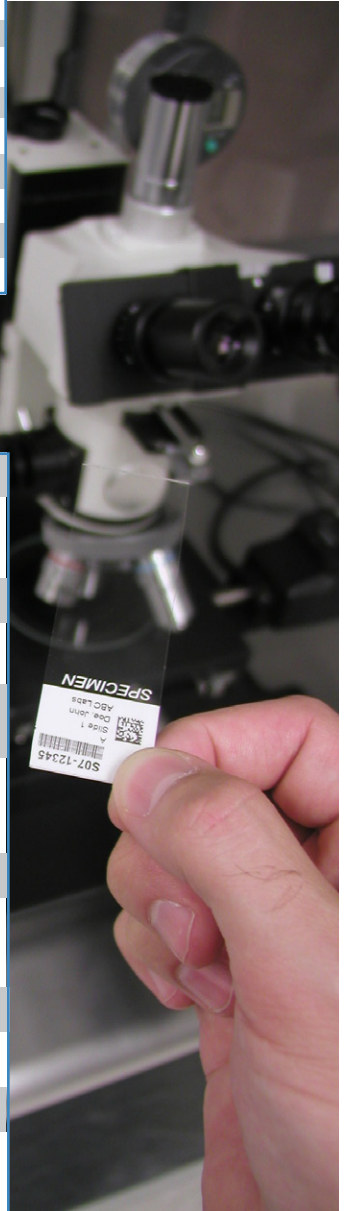

**pslim**  
Laboratory Slide Printer

**ACCUp**lace

(954) 791-1500  
1800 NW 69th Avenue  
Plantation, FL 33313  
accuplace.com